

Ethylene-Oxide

Methanol

Ethanol

Acetone


IPA

Acetaldehyde

1,3-Butadiène

airmoVOC C1C3

Analysis of light oxygenated VOC's in air or pure gas

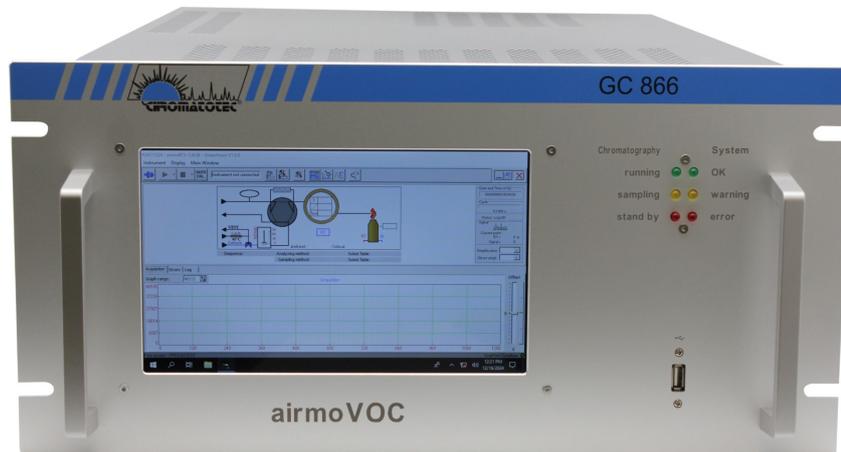

Model: A13000 (5U) computer in option

Environment:

- Urban/Non urban area pollution control
- Indoor measurements
- PAMS /CE analysis
- Out of plant

Process:

- Industrial Hygiene
- Fence line

Toxicity:

- VOC screening in food and beverage industry

SCAN or CLICK ME



Chromatotec® is specialized in VOC, Sulfur and permanent gases analysis at trace and ultra trace levels (ppm, ppb, ppt). Please visit our website for more details.

Updated: September 2025

airmoOVOC C1C3

Analysis of light oxygenated VOC's in air or pure gas



Principle:

This instrument uses a 6-port valve with a single sorbent trap, and metallic capillary columns.

It has a minimum detection level of as low as 0.1 ppb in ambient air. The airmo OVOC C1C3 uses a 3 phases trap for the analysis of light oxygenated VOCs.

Miniaturization, sensitivity, mobility and flexibility are its main features. Everything from the sample port up to the data storage is integrated in a 19" - rack 5U.

Programmed gradient temperature of the oven and pressure/flow control of the carrier gas.

Before final delivery the analyser is tested for three weeks by the quality control department.

The Vistachrom software enables the user to view and store data on an industrial computer. It provides comfortable user friendly utilities to recalculate, calibrate and export data and to configure the measurement.

The software allows the calculation of retention time, area, mass or concentration profiles.

Calib is a very important option for good result in ppb range, H2 is necessary.

We use permeation tubes for calibration

Product technical specifications:

Light VOC analysis:

- Methanol/Ethanol
- Ethylene-Oxide
- acetone
- IPA
- Acetaldehyde
- 1,3-Butadiène

Detection limit:

- Less than 0.1 ppb

Detection range:

- 0 - 100 ppb / amplification 3 in standard
- 0 - 1000 ppb
- 0 - 10 000 ppb

Relative Standard Deviation:

- < 0.3 % over 48h (RT)
- < 2 % over 48 h (Conc.)

Results:

- Data storage
- 4-20mA current output (option)
- MODBUS / JBUS or MGS1 communication protocol (option)

Cycle time:

- 30 min as standard

Gas supply:

- H2 (FID): 30 ml/min
- Air (FID): 180 ml/min
- Carrier gas H2
- Air (FID): 180 ml/min (inlet 3 bars ; 1/8")
- Sample inlet (vacuum pump) 1/4"
- Pneumatic valve : 90ml/commutation

Sample volume:

- 20 to 240 ml (programmable) in standard

Power supply:

- Main: (230V / 115V 50 Hz/60Hz)
- Battery: 24V (option)

Electrical consumption:

- Mean: 180 VA, peak 400 VA

Dimensions and weight:

- Rack: 19" - (5U)
- Height: 222mm
- Width: 482mm
- Depth: 660mm
- Net weight: 20 Kg

[CLICK HERE FOR ADDITIONAL DIMENSIONS DETAILS](#)

Options:

- MODBUS / JBUS or MGS1
- Analog output 4-20 mA or 0-10 V or alarms
- Automatic validation and calibration
- 24 V power supply
- Hydrogen and zero air generators for transportable analyzers
- Multiple Stream selector (2 to 6)
- Nitrogen generator

To order:

airmo OVOC C1C3 (5U)
Needs a supervisor
Supervisor (5U)

Model:

A13000
XXX022

Chromatotec® is continuously improving its products, therefore these specifications are subject to change without notice

To contact us: sales@chromatotec.com

NORTH AMERICA
Houston - USA

EUROPE
Bordeaux - FRANCE

ASIA
Beijing - CHINA