

Zero VOC air

Flow security

Tank

Air generator

Air dryer

VOC catalyser

# airmoPURE D 45 PSI

Zero air generator



Model: XXX031-D

Zero VOC air production for the analyzer or for carrier gas

*Chromatotec® is specialised in VOC, Sulfur and permanent gases analysis at trace and ultra trace levels (ppm, ppb, ppt).  
Please visit our website for more details.*

Updated: March 2024

# airmoPURE D 45 PSI

Zero air generator



## The airmoPure provides zero air for:

- FID detector: 180 ml/min
- Pneumatic valve: 60 ml by commutation
- Permapure: 150 ml/min for C2C6: **Total flow for airmoOzone = 580 ml/min**
- Cooling box (Peltier effect) for C2C6: 13 ml/min
- Dilution of permeation tube for calibration : 50 ml/min
- Carrier gas for MEDOR®

## The airmoPure includes:

- A compressor station fixed with a **10-liter reservoir** at 7 bars which delivers outlet pressure at 3 or 4 bars with **protection**
- A membrane air dryer with automatic purge. The air is dried before being compressed in the cylinder.
- A heated catalyser to eliminate hydrocarbons at ppt level
- Connections between the airmoPure and the instrument

## Example with 3 bars inlet:

- 180 ml/min for FID
- 50 ml/min in continuous for Calibration
- Dilution: 230 ml/min in discontinuous mode (during calibration only)
- Benzene permeation tube (32 ng/min  $\pm$  10% at 45°C) diluted at 230 ml/min = 139  $\mu$ g/m<sup>3</sup> or 37 ppb
- Zero air for analysis in continuous mode (airmoBTX /CALIB): 180 + 50  $\approx$  230 ml/min
- Zero air for analysis in discontinuous (airmoBTX/CALIB): 180 + 230  $\approx$  410 ml/min

## Product technical specifications:

### Zero air total hydrocarbon concentration:

- < 0.1 ppb NMTHC in ambient air without pollution

### Dew point:

- < -15°C°

### Zero air maximum flow:

- 1 l/min in discontinuous mode for airmoOzone

### Noise:

- 55 dB(A)

### Zero air start up time:

- 45 min

### Power supply:


- 230V / 110V / 115V and 50 Hz/60Hz

### Electrical consumption:

- 310 VA

### Dimensions and weight:

- Height: 370 mm
- Width: 370 mm
- Depth: 370 mm
- Net Weight: 17 kg

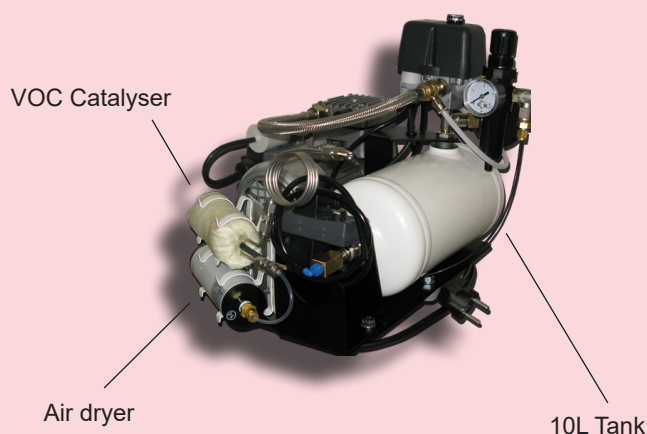

**To order:**  
airmoPURE D 45 PSI

**Model:**  
XXX031-D

Chromatotec® is continuously improving its products, therefore these specifications are subject to change without notice



To contact us: [sales@chromatotec.com](mailto:sales@chromatotec.com)

**NORTH AMERICA**  
Houston - USA

**EUROPE**  
Bordeaux - FRANCE

**ASIA**  
Beijing - CHINA

**WWW.CHROMATOTEC.COM**